



First Committee Draft (CD) – Marked version

Project: New Recommendation
Title: Ophthalmic instruments – Non-contact tonometers
Part 1: Technical and metrological requirements
Date: 10 August 2020
Document number: TC18_P2_N008
Supersedes document: -
Project Group: OIML TC 18/p 2
Convenership: Germany
Convener: Dr Dana Maria Rosu

Circulated to P- and O-members and liaison international bodies and external organizations for:

☐ Discussion at (date and place of meeting):

☒ Comments by: 10 November 2020

☐ Vote (P-members only) and comments by:

Contents (to be finalised later)

Forward	4
1 Introduction	5
2 Scope	5
3 Terms and definitions	5
3.1 intraocular pressure (IOP)	6
3.2 non-contact tonometer (NCT)	6
3.3 reference tonometer	6
3.5 Hand-held	6
4 Description of the category of instrument	6
5 Units of measurement	6
6 Metrological requirements	7
6.1 IOP measuring range and resolution	7
6.2 Environmental conditions	7
6.3 Clinical accuracy	7
6.4 Stability	8
7 Technical requirements	9
7.1 Mechanical strength of non-contact tonometers	9
8 Metrological controls	9
8.1 Inscription of the device	9
8.2 Accompanying documents	9
8.3 Verification	10
Mandatory Annex – IOP testing equipment	11
A.1 Requirements for IOP testing equipment	11
A.1.1 IOP testing equipment in use	11
A1.2 Storage of IOP testing equipment	11
A1.3 Mechanical stability	12

A1.4 Mechanical strength of IOP testing equipment	12
Informative Annex	13

Forward

The International Organization of Legal Metrology (OIML) is a worldwide, intergovernmental organization whose primary aim is to harmonize the regulations and metrological controls applied by the national metrological services, or related organizations, of its Member States.

The main categories of OIML publications are:

- **International Recommendations (OIML R)**, which are model regulations that establish the metrological characteristics required of certain measuring instruments and which specify methods and equipment for checking their conformity. OIML Member States shall implement these Recommendations to the greatest possible extent;
- **International Documents (OIML D)**, which are informative in nature and which are intended to harmonize and improve work in the field of legal metrology;
- **International Guides (OIML G)**, which are also informative in nature and which are intended to give guidelines for the application of certain requirements to legal metrology; and
- **International Basic Publications (OIML B)**, which define the operating rules of the various OIML structures and systems.

OIML Draft Recommendations, Documents and Guides are developed by Project Groups linked to Technical Committees or Subcommittees which comprise representatives from the Member States. Certain international and regional institutions also participate on a consultation basis. Cooperative agreements have been established between the OIML and certain institutions, such as ISO and the IEC, with the objective of avoiding contradictory requirements. Consequently, manufacturers and users of measuring instruments, test laboratories, etc. may simultaneously apply OIML publications and those of other institutions.

International Recommendations, Documents, Guides and Basic Publications are published in English (E) and translated into French (F) and are subject to periodic revision.

Additionally, the OIML publishes or participates in the publication of **Vocabularies (OIML V)** and periodically commissions legal metrology experts to write **Expert Reports (OIML E)**. Expert Reports are intended to provide information and advice, and are written solely from the viewpoint of their author, without the involvement of a Technical Committee or Subcommittee, nor that of the CIML. Thus, they do not necessarily represent the views of the OIML.

The present document is the first Committee Draft (1CD was developed by Project Group 2 of OIML TC 18 Medical measuring instruments.) and was drawn up on the basis of the conclusions of comments from member nations on the Working Draft circulated on the 15th of January 2020. This publication was developed by Project Group 2 of OIML TC 18 Medical measuring instruments.

OIML Publications may be downloaded from the OIML website in the form of PDF files. Additional information on OIML Publications may be obtained from the Organization's headquarters:

Bureau International de Métrologie Légale
11, rue Turgot - 75009 Paris – France
Telephone: 33 (0)1 48 78 12 82
Fax: 33 (0)1 42 82 17 27
E-mail: biml@oiml.org
Internet: www.oiml.org

Ophthalmic instruments – Non-contact tonometers

Part 1 – Metrological and technical requirements

1 Introduction

Non-contact tonometers have been used for many years to measure the intraocular pressure. Accurate measurement of the IOP is important as chronically high IOP is considered the main risk factor for glaucoma development and progression.

The Recommendation references the international standard ISO 8612 Ophthalmic instruments – Tonometers with respect to the clinical investigation of non-contact tonometers.

Note: ISO 8612 uses the term “design compliance testing” for clinical investigation.

To ensure the traceability of non-contact tonometers, this Recommendation proposes a series of requirements and test procedures meant to compare the performance of individual non-contact tonometers to the, in its performance and safety components, identically manufactured clinically tested device in conformity with ISO 8612 (commonly known as golden unit).

Additional requirements concerning environmental influences on the performance and mechanical strength are specified.

2 Scope

This Recommendation specifies minimum requirements and test procedures for the design compliance, construction and performance of non-contact tonometers (NCT), which are used for the estimation of the intraocular pressure (IOP) of human subjects in clinical applications.

The requirements and test procedures for applanation and impression tonometers are described in the OIML Recommendation R 145.

Note: Non-contact tonometers do not measure directly the intraocular pressure; the IOP is estimated by investigating the response of the eye to an air pulse generated by a non-contact tonometer.

In the mandatory [Annex B](#) of this Recommendation, minimum requirements and test procedures for the design compliance, construction and performance of IOP testing equipment are specified.

3 Terms and definitions

For the purposes of this Recommendation the following terms and definitions apply.

3.1 intraocular pressure (IOP)

pressure within the eye front chamber; units given in millimetres of mercury (mmHg) or kilopascals (kPa)

3.2 non-contact tonometer (NCT)

tonometer that uses an air pulse to determine the IOP, without mechanical contact with the eye.

3.3 reference tonometer

~~A~~ Goldman applanation tonometer that complies with Annex A of ISO 8612.

3.4 IOP testing equipment

~~E~~quipment able to simulate intraocular pressure and which is used for the testing and verification of non-contact tonometers.

Note: The most common types of IOP testing equipment for non-contact tonometers are described in the ~~Informative Annex~~ informative Annex B.

3.5 Hand-held

~~refers to a~~ device intended to be supported by the hand during normal use.

4 Description of the category of instrument

Non-contact tonometers use the applanation principle. The applanation of a predefined area of the cornea is achieved by the application of an air pulse. An optical sensor registers when the applanation is achieved and the IOP is calculated based on the force needed to applanate the cornea.

The main components of a non-contact tonometer (air-puff tonometer) are: a pneumatic system that generates a rapid puff of air directed to the eye surface, a light emitter that sends for instance a beam of light under a specified angle and a light detector collecting the light reflected by the cornea at a specified angle of incidence.

5 Units of measurement

The IOP shall be indicated in kilopascals (kPa) or in millimetres of mercury (mmHg).

6 Metrological requirements

6.1 IOP measuring range and resolution

The non-contact tonometers shall be capable of indicating IOP at least in the range 0.0 kPa ~~to~~ 5.3 kPa (0 mmHg ~~to~~ 40 mmHg). If the manufacturer claims a broader range, the compliance with this ~~R~~Recommendation shall be proven for the entire range.

The resolution of the display shall be at least 0.1 kPa (1 mmHg).

Check compliance with the requirements consulting the instructions for use or by confirmation by the manufacturer.

6.2 Environmental conditions

6.2.1 Devices in use

Non-contact tonometers shall comply with all the requirements specified in 6.4.1 at

- temperatures between 10 °C and 40 °C,
- relative humidity between 10 % and 90 % (non-condensing).

Testing shall be carried out in accordance with Clause 1 of R XXX-2.

6.2.2 Storage of non-contact tonometers

Non-contact tonometers shall comply with all requirements specified in 6.4.1 after storage at

- temperatures between –20 °C and 60 °C,
- relative humidity between 10 % and 90 % (non-condensing).

Testing shall be carried out in accordance with Clause 2 of R XXX-2.

6.3 Clinical investigation

The manufacturer shall demonstrate, on the basis of design compliance testing as specified in Clause 4.2 of ISO 8612, that at least 95 % of the paired differences between the reference tonometer reading and the test tonometer reading, meet the requirements given in Table 1.

Table 1
The maximum permissible measurement error
for the paired differences reference tonometer - test tonometer

IOP range			maximum permissible measurement error	
	kPa	mmHg	kPa	mmHg
low	0.93 to 2.13	7.0 to 16.0	± 0.67	± 5.0
intermediate	2.13 to < 3.07	> 16.0 to < 23.0	± 0.67	± 5.0
high	≥ 3.07	≥ 23.0	± 0.67	± 5.0

Note: ISO 8612 uses the term “tolerance” as “maximum permissible measurement error”. This recommendation opts for the use of the term “maximum permissible measurement error” as recommended in the International vocabulary of metrology - Basic and general concepts and associated terms (VIM, 4.26).

6.4 Stability

The repeatability test in Clause 6.4.1 shall be performed as the initial test (prior to all other tests required by this Recommendation).

The intermediate measurement precision test in Clause 6.4.2 shall be performed as the final test (following all other tests required by this Recommendation).

6.4.1 Repeatability

~~40~~**Ten** repeated measurements shall be performed for one value in each of the three IOP ranges specified in the Table 1. This value is specified by the manufacturer for the IOP testing equipment used.

For any temperature-humidity combination in the temperature range from 10° C to 40 °C and the relative humidity range from 10 % to 90 % (non-condensing), the maximum permissible measurement error Δx (the difference between the indicated value of the device under test and the reference value) of the arithmetic mean and the standard deviation $s(\Delta)$ shall comply with the following requirements:

- for the reference value in the low range: $|\Delta x_{low}| \leq 0.13 \text{ kPa (1.0 mmHg)}$; $s(\Delta x_{low}) \leq 0.13 \text{ kPa (1.0 mmHg)}$
- for the reference value in the intermediate range: $|\Delta x_{intermediate}| \leq 0.20 \text{ kPa (1.5 mmHg)}$; $s(\Delta x_{intermediate}) \leq 0.20 \text{ kPa (1.5 mmHg)}$
- for the reference value in the high range: $|\Delta x_{high}| \leq 0.27 \text{ kPa (2.0 mmHg)}$; $s(\Delta x_{high}) \leq 0.27 \text{ kPa (2.0 mmHg)}$

•

Testing shall be carried out in accordance with Clause 1 of R XXX-2

Note: This Clause does not require new measurements if the measurements from Clause 6.2 have already been recorded.

6.4.2 Intermediate measurement precision

The repeatability test in Clause 6.4.1 shall be repeated and complied with.

Testing shall be carried out in accordance with Clause 3 of R XXX-2

7 Technical requirements

7.1 Mechanical strength of non-contact tonometers

Hand-held non-contact tonometers shall operate correctly following a free fall from a height of 1.0 ± 0.1 m as specified in Clause 4 of R xxx-2.

~~Non hand-held non contact tonometers shall operate correctly following a free fall from a height of 5 cm as specified in Clause 4 of Rxxx-2.~~

~~In all cases, e~~Correct operation following the falls shall be verified by checking that the tonometer still complies with the requirements in Clause 6.4.

Testing shall be carried out in accordance with Clause 4 of R XXX-2

8 Metrological controls

8.1 Inscription of the device

Each tonometer shall bear the following information:

- name of manufacturer or ~~trade mark~~trademark;
- serial number;
- measuring unit;
- type approval number (if applicable).

Testing shall be carried out by visual inspection.

8.2 Accompanying documents

The manufacturers shall provide a user's instruction manual, which includes:

- the description of the device and its accessories.
- instructions for use and maintenance.
- information concerning clinical investigation according to ISO 8612 Ophthalmic instruments – Tonometers.
- reference to the current Recommendation.

Testing shall be carried out by visual inspection.

8.3 Verification

In the case of non-contact tonometers, traceability to national standards is technically not feasible. Instead, traceability is ensured via a transfer standard, to an identically manufactured device in its performance and components and clinically tested in conformity with ISO 8612. IOP testing equipment allowing the simulation of IOP measurements is used as transfer standards in the (re)verification of non-contact tonometers.

The verifications shall be periodically performed and after repair. The reverification period shall not exceed 2two years.

The (re)verification shall be carried out in accordance to Clause 1 of R XXX-2 at any room temperature between 10 °C -and 40 °C and relative humidity between 10 % -and 90 %.

Note: Examples for such IOP testing equipment can be found in Annex 4B. The use of different IOP testing equipment is permissible as long as the manufacturer can prove compliance with the current Recommendation.

~~Mandatory Annex A-~~

IOP testing equipment

(Mandatory)

A.1 Requirements for IOP testing equipment

The IOP testing equipment shall prove sufficiently small error limits and measurement uncertainties. The error limits are considered to be sufficiently small if they do not exceed one third of the error limits of the NCT to be tested. Manufacturers may exceed these requirements, i.e. to specify smaller permissible limits for their IOP testing equipment. Manufacturer specifications below the state of the art shall not be acceptable for IOP testing equipment intended for metrological control.

Note: Examples for such IOP testing equipment can be found in Annex ~~1B~~. The use of different IOP testing equipment is permissible as long as the manufacturer can prove compliance with the current Recommendation.

The IOP testing equipment shall be periodically verified according to A.1.1. The reverification period is ~~1~~one years unless otherwise specified by the manufacturer, but not more than 3three years.

A.1.1 IOP testing equipment in use

The IOP testing equipment shall comply with all the requirements specified in this Recommendation at

- temperatures between 10 °C and 40 °C,
- relative humidity between 10 % and 90 % (non-condensing).

Check compliance with Clause 6.4.1.

Testing shall be carried out in accordance with Clause 1 of R XXX-2

A.1.2 Storage of IOP testing equipment

The test equipment shall comply with all the requirements specified in this Recommendation at

- temperatures between 10 °C and 40 °C,
- relative humidity between 10 % and 90 % (non-condensing).

Check compliance with Clause 6.4.1.

Testing shall be carried out in accordance with Clause 2 of R XXX-2

A1.3 Mechanical stability

Perform at least 1 000 repeated measurements with the IOP testing equipment, equally distributed in the three measurement ranges specified in Table 1.

Check compliance with 6.4.1.

Testing shall be carried out in accordance with Clause 4 of R XXX-2

A1.4 Mechanical strength of IOP testing equipment

The [hand-held](#) IOP testing equipment shall operate correctly following a free fall from a height of $1.0 \text{ m} \pm 0.1 \text{ m}$ onto a $50 \text{ mm} \pm 5 \text{ mm}$ thick hardwood board (density $> 600 \text{ kg/m}^3$) lying flat on a concrete or similar rigid base. This requirement shall apply to a fall from a starting position in any of the three axes of orientation.

IOP testing equipment that cannot comply with these requirements shall be marked accordingly and shall be subject to performance checks prior to each use.

Check compliance with 6.4.1.

Testing shall be carried out in accordance with Clause 4 of R XXX-2

~~Informative~~ Annex B

IOP testing equipment to simulate the measurement of the intraocular pressure

(Informative)

The IOP testing equipment to simulate the measurement of the intraocular pressure shall fulfil at least the following requirements:

- The non-contact tonometer to be tested in the same configuration as it is used on human subjects.
- The IOP testing equipment to be stable (see Clause 10).

Note: The procedure on how to attach the IOP testing equipment to the non-contact tonometer should be part of the testing.


Those components important for the simulation to be traceable to national standards.

Different set-ups to test non-contact tonometers exist and new ones are being developed. As a presentation of all available IOP testing equipment is not feasible, this ~~¶~~Recommendation will concentrate on the three most commonly used set-ups:

- PTB flapper
- mechanical phantom eye (“rubber eye”)
- electronic phantom eye

PTB flapper

The German Metrology Institute PTB (Physikalisch-Technische Bundesanstalt) has developed in the 1980s a mechanical device able to simulate an IOP measurement. The scheme of this IOP testing equipment is shown in Figure 1. This device simulates the applanation of a human eye by means of the kinematic movement of a small circular mirror (7). The device uses a mirror and lever system, wherein the tonometer air pulse is directed at the mirror mounted on a lever to angularly displace the lever about a pivot axis. By changing the torque, different IOP values can be simulated.

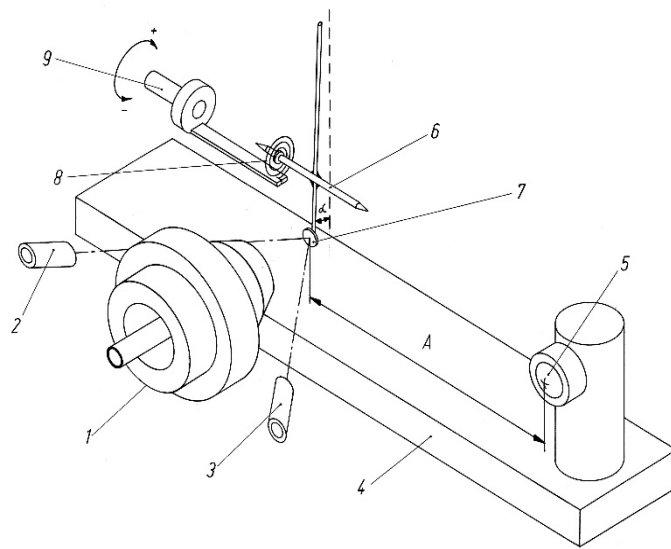

Where:

- 1: nozzle of the non-contact tonometer
- 2, 3: light source and light detector of the NCT
- 4: mount
- 5: lens to adjust the distance non-contact tonometer – mirror
- 6: rotating axis
- 7: mirror
- 8: spring
- 9: adjustment of the torque
- A: distance to move the test jig distant adjustment to the dynamic test position

Figure 1 - Schematic drawing of the PTB flapper

Note: The PTB flapper requires special mounting to be attached to the non-contact tonometer.

Mechanical phantom eye (“rubber eye”)

The mechanical phantom eye shall simulate the shape and the biomechanical properties of the cornea, it is made of soft silicon rubber and fabricated with different thicknesses. The aging effect of the elastic material shall be investigated, and periodic checks are strongly recommended. An example of a mechanical phantom eye is shown in Figure 2.

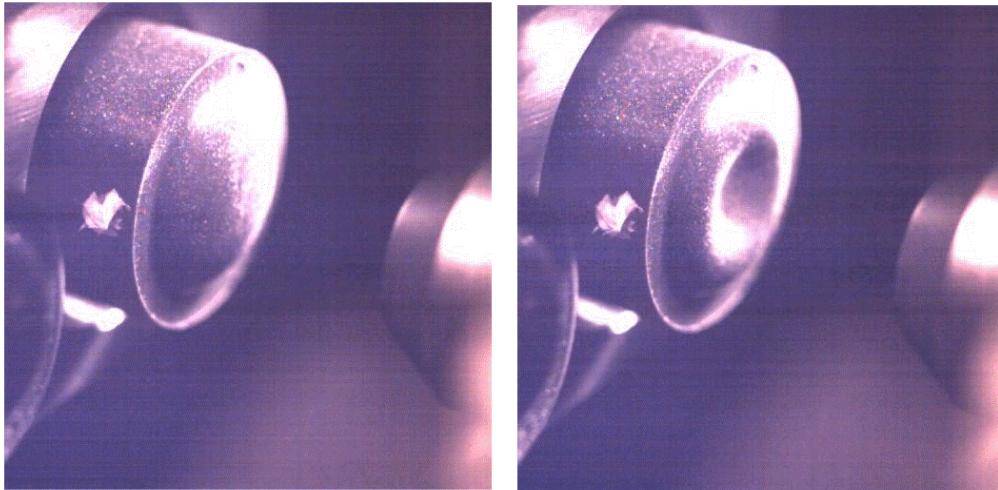



Figure 2 - Representation of a mechanical phantom eye
(in both cases shown on the left side of the picture)

left: in the resting phase; right: during the measurement

Note: Usually this mechanical phantom eye can be attached to the non-contact tonometer easily.

Electronic phantom eye

The electronic phantom eye (Figure 3) is a device, which is made from rigid material. Its exterior design resembles a cornea. Inside it has a pressure transducer and a light source. Depending on the pressure applied by the air puff, measured by the pressure transducer of the electronic phantom eye, the light source flashes. This flash simulates the reflection signal expected by the non-contact tonometer. The electronic phantom eye is mechanically stable and easy to attach to the non-contact tonometer.

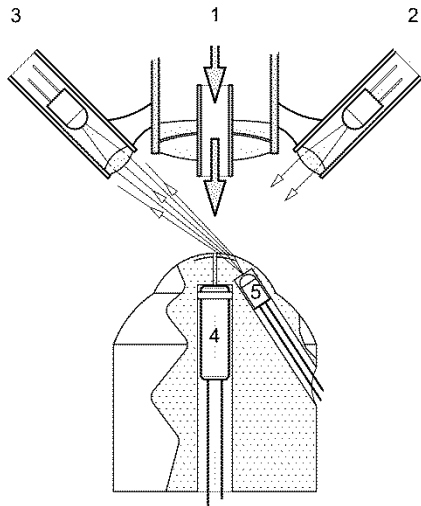

Where

- 1: nozzle of the non-contact tonometer
- 2, 3: light source and light detector of the NCT
- 4: pressure transducer of the electronic phantom eye
- 5: light source of the electronic phantom eye

Figure 3 - Schematic drawing of an electronic phantom eye